

10/531844
Rec'd PCT/PTO 18 APR 2005
PCT/2003/004542

#2



INVESTOR IN PEOPLE

The Patent Office
Concept House
Cardiff Road
Newport
South Wales
NP10 8QQ

REC'D 04 DEC 2003

WIPO PCT

**PRIORITY
DOCUMENT**
SUBMITTED OR TRANSMITTED IN
COMPLIANCE WITH RULE 17.1(a) OR (b)

I, the undersigned, being an officer duly authorised in accordance with Section 74(1) and (4) of the Deregulation & Contracting Out Act 1994, to sign and issue certificates on behalf of the Comptroller-General, hereby certify that annexed hereto is a true copy of the documents as originally filed in connection with the patent application identified therein.

In accordance with the Patents (Companies Re-registration) Rules 1982, if a company named in this certificate and any accompanying documents has re-registered under the Companies Act 1980 with the same name as that with which it was registered immediately before re-registration save for the substitution as, or inclusion as, the last part of the name of the words "public limited company" or their equivalents in Welsh, references to the name of the company in this certificate and any accompanying documents shall be treated as references to the name with which it is so re-registered.

In accordance with the rules, the words "public limited company" may be replaced by p.l.c., plc, P.L.C. or PLC.

Re-registration under the Companies Act does not constitute a new legal entity but merely subjects the company to certain additional company law rules.

Signed

R. Mahoney

Dated 10 November 2003

An Executive Agency of the Department of Trade and Industry

BEST AVAILABLE COPY

Patents Form 1/77

Patents Act 1977
(Rule 16)

The
Patent
Office

THE PATENT OFFICE
H

26 OCT 2002

NEWPORT

2800002 E 58983-1/002934
P01/7700 0.00-0224997.7

Request for grant of a patent

(See the notes on the back of this form. You can also get
an explanatory leaflet from the Patent Office to help
you fill in this form)

The Patent Office

Cardiff Road
Newport
Gwent NP9 1RH

1. Your reference 100874

2. Patent application number 0224997.7
(The Patent Office will fill in this part)

26 OCT 2002

3. Full name, address and postcode of the or of
each applicant (underline all surnames) AstraZeneca AB
S-151 85 Sodertalje
Sweden

Patents ADP number (if you know it)

If the applicant is a corporate body, give the
country/state of its incorporation Sweden

4. Title of the invention METHOD

5. Name of your agent (if you have one) Neil Godfrey Alasdair Phillips

"Address for service" in the United Kingdom
to which all correspondence should be sent
(including the postcode)

AstraZeneca UK Limited
Global Intellectual Property
Mereside, Alderley Park
Macclesfield
Cheshire SK10 4TG

Patents ADP number (if you know it)

| 6. If you are declaring priority from one or more earlier patent applications, give the country and the date of filing of the or of each of these earlier applications and (if you know it) the or each application number | Country | Priority application number (if you know it) | Date of filing (day / month / year) |
|--|---------|---|--|
| | | | |

| 7. If this application is divided or otherwise derived from an earlier UK application, give the number and the filing date of the earlier application | Number of earlier application | Date of filing (day / month / year) |
|--|-------------------------------|--|
| | | |

8. Is a statement of inventorship and of right
to grant of a patent required in support of
this request? (Answer 'Yes' if:
a) any applicant named in part 3 is not an inventor, or
b) there is an inventor who is not named as an
applicant, or
c) any named applicant is a corporate body.
See note (d))

Patents Form 1/77

Patents Form 1/77

9. Enter the number of sheets for any of the following items you are filing with this form. Do not count copies of the same document

Continuation sheets of this form

Description

7 ✓

Claim(s)

1 ✓

Abstract

Drawing(s)

8 + 8

10. If you are also filing any of the following, state how many against each item.

Priority documents

Translations of priority documents

Statement of inventorship and right to grant of a patent (Patents Form 7/77)

Request for preliminary examination and search (Patents Form 9/77)

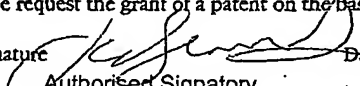
Request for substantive examination (Patents Form 10/77)

Any other documents (please specify)

11.

I/We request the grant of a patent on the basis of this application.

Signature



Date

Authorised Signatory

25/10/2002

12. Name and daytime telephone number of person to contact in the United Kingdom

Jennifer C Bennett - 01625 230148

Warning

After an application for a patent has been filed, the Comptroller of the Patent Office will consider whether publication or communication of the invention should be prohibited or restricted under Section 22 of the Patents Act 1977. You will be informed if it is necessary to prohibit or restrict your invention in this way. Furthermore, if you live in the United Kingdom, Section 23 of the Patents Act 1977 stops you from applying for a patent abroad without first getting written permission from the Patent Office unless an application has been filed at least 6 weeks beforehand in the United Kingdom for a patent for the same invention and either no direction prohibiting publication or communication has been given, or any such direction has been revoked.

Notes

- If you need help to fill in this form or you have any questions, please contact the Patent Office on 0645 500505.
- Write your answers in capital letters using black ink or you may type them.
- If there is not enough space for all the relevant details on any part of this form, please continue on a separate sheet of paper and write "see continuation sheet" in the relevant part(s). Any continuation sheet should be attached to this form.
- If you have answered 'Yes' Patents Form 7/77 will need to be filed.
- Once you have filled in the form you must remember to sign and date it.
- For details of the fee and ways to pay please contact the Patent Office.

METHOD

The present invention relates to improved screening assays and in particular to the use of activator-independent forms of the murein biosynthesis enzyme Mur D, such as from
5 *Enterococcus faecalis* (*E. faecalis*). Such screening assays are used to identify and characterize modulators of the Mur D enzyme.

Interest has been shown in the murein (Mur) biosynthesis pathway in bacteria, this is a key component of bacterial cell wall synthesis. Enzymes in this pathway are potential targets for broad-spectrum and selective antibacterial agents.

10 Mur D (UDP-N-acetylmuramyl-L-alanine:D-glutamate ligase) is the bacterial D-glutamate-adding enzyme which catalyses the attachment of D-glutamate to a cytoplasmic peptidoglycan precursor, UDP-N-acetylmuramyl-L-alanine. This reaction results in the formation of a peptide linkage between the amino function of D-glutamate and the carboxyl terminus of UDP-N-acetylmuramyl-L-alanine. A stoichiometric consumption of ATP
15 supplies the energy needed for this peptide bond formation resulting in generation of ADP and orthophosphate.

Walsh et al (Journal of Bacteriology, Sept 1999, 181, No.17, 5395-5401) have examined the biochemical properties of the Mur D enzyme from two gram-negative bacteria ie. *Escherichia coli*, and *Haemophilus influenzae*, and two gram-positive bacteria ie.

20 *Enterococcus faecalis* and *Staphylococcus aureus*. They established data regarding the biochemical properties of these enzymes and discussed similarities and differences between them, in particular salt-activation of the gram-negative bacteria. They report that the differences observed between the gram-positive and the gram-negative bacteria indicate that the two gram-negative bacteria may apply a more stringent regulation of cell wall
25 biosynthesis at the early stage of the peptidoglycan biosynthesis pathway than do the two gram-positive bacteria. However, the skilled reader is aware that the substrate purification procedure used by Walsh et al cannot remove all salts that may function as enzyme activators. Therefore it is not possible to draw meaningful conclusions as to the salt-dependency or otherwise of the two gram-positive bacteria.

30 Indeed, we have now found that the Mur D enzyme from the gram-positive bacterium *Staphylococcus aureus* is also salt-activated. This led to our discovery that the *E. faecalis* Mur D enzyme has unique properties which make it possible to devise improved screening assays using an activator-independent Mur D enzyme.

Therefore in a first aspect of the present invention, we provide the use of an activator-independent MurD enzyme in a screening assay to identify inhibitors of the enzyme, which assay comprises contacting the enzyme with a test compound in the presence of an enzyme substrate and appropriate buffers etc. and detecting any modulation of enzyme activity by the test compound.

By "activator-independent" we mean that the Mur D enzyme is not activated by salt species normally associated with the substrate (here D-glutamic acid) or other assay components. Monovalent cations such as ammonium (NH₄⁺) and potassium (K⁺) are particular salt species that activate Mur D. We note that if different amounts of substrate are used in an assay e.g. for K_m determinations or mode of inhibition studies, the amount of the activating ammonium ions is not constant. It can therefore not be clearly identified if an activity increase is due to activation or due to an increase of substrate concentration. Contrary to the results reported by Walsh et al we have found that the *Staphylococcus aureus* Mur D enzyme is activated by NH₄⁺ and K⁺ cations (cf. Figure 2).

Our analysis shows that activator-dependent forms of the MurD enzyme have the following common amino acid residues ie. G96, A112, A116, V126, L129, M133, G296, P298 and V422. The indicated positions are based on the *E. faecalis* Mur D sequence as set out in Figure 8 (and corresponding sequence alignments).

Therefore in a further aspect of the invention we provide the use of an activator-independent MurD enzyme which comprises a Mur D amino acid sequence wherein one or more of the amino acid residues at the positions given above is not as indicated for that particular amino acid. More conveniently, at least two, or at least three, or at least four, or at least five, or at least six, or at least seven, or at least 8, or all nine of the amino acid residues are not as indicated.

The activator-independent Mur D may conveniently comprise an amino acid sequence comprising one of more of the following residues ie. K96, C112, G116, T126, M129, L133, N296, S298 and I422. . More conveniently, at least two, or at least three, or at least four, or at least five, or at least six, or at least seven, or at least 8, or all nine of the amino acid residues are as indicated.

The activator-independent MurD enzyme conveniently comprises the following *E faecalis* amino acid sequence

MKKITTYQNK KVLVLGLAKS GVSAKLLHE LGALVTVNDA KQFDQNPDAQ
DLLTLGIRVV TGGHPIELLD EEFELIVKNP GIPYTNPLVA EALTRKIPII

TEVELAGQIA ECPVIGITGT NGKTTTTTMI GLLLNADRTA GEARLAGNIG
 FPASTVAQEA TAKDDLVMEL SSFQLMGIET FHPQIAVITN IFEAHLDYHG
 SRKEYVAAKW AIQKNMTVED TLILNWNQVE LQTLAKTTAA NVLPFSTKEA
 VEGAYLLDGK LYFNEEYIMP ADELGIPGSH NIENALAAIC VAKLKNVSNV
 5 QIRQTLKNFS GVPHRTQFVG EVQRRFYND SKATNILATE MALSGFDNQK
 LLLLAGGLDR GNSFDELVPA LLGLKAIVLF GETKEKLAEA AKKANIELIL
 FAENVQTAVT IAFDYSEKDD TILLSPACAS WDQYPNFEVR GEAFMQAVQQ
 LKESEM

or a Mur D amino acid sequence having at least 85% homology, such as 90 or 95%

10 homology, therewith.

Alternatively the activator-independent Mur D amino acid sequence has at least 70% sequence identity with the above amino acid sequence.

Whilst we do not wish to be limited by theoretical considerations we believe that the activator-independent Mur D enzyme may have up to 16, such as up to 12, 10, 8, 6, 4, or 2,
 15 amino acids removed from the N-terminus and/or up to 12, such as up to 10, 8, 6, 4, or 2, amino acids removed from the C-terminus of the enzyme.

Activator-independent Mur D sequences that do not correspond to published Mur D enzyme sequences are novel and represent a further aspect of the present invention.

Any convenient screening assay format may be used. By way of non-limiting
 20 example we disclose the following:

The enzyme substrates are conveniently UDP-MurNac-L-Ala, D-Glutamate and ATP.

The enzyme is conveniently pre-incubated with the test compound to allow inhibitors to bind to the enzyme. This may allow the detection of inhibitors with a slow binding mode to the enzyme or allow detection of specifically modifying inhibitors that may be out-competed
 25 by the substrates.

Any appropriate buffer can be used that has a pKa in the active pH range of *E. faecalis* MurD (pH 7.0 – 10.0) Examples of convenient buffers include buffers that do not contain phosphate such as GOOD Buffer i.e. Tris or Hepes (Good, et al. (1966) Biochemistry, 5, 467-477). Detection of modulation of enzyme activity may be effected using any convenient
 30 detection system, such as those which include a colour change eg. using malachite green. These include absorbance spectrophotometers, absorbance plate reader or any other instrument that can determine the absorption of a solution between for example 400 and 800 nm

Modulation of enzyme activity can be enzyme inhibition or activation or enhancement of enzyme activity, conveniently enzyme inhibition.

Appropriate control reactions are conveniently performed to determine if a chemical compound interferes with the detection system and/or has an absorbance at the detection

5 wavelength.

The test compound is any convenient compound that may be useful in pharmaceutical research. It may be a polypeptide of equal to or greater than, 2 amino acids such as up to 6 amino acids, up to 10 or 12 amino acids, up to 20 amino acids or greater than 20 amino acids such as up to 50 amino acids. For drug screening purposes, preferred compounds are

10 chemical compounds of low molecular weight and potential therapeutic agents. They are for example of less than about 1000 Daltons, such as less than 800, 600 or 400 Daltons in weight. If desired the test compound may be a member of a chemical library. This may comprise any convenient number of individual members, for example tens to hundreds to thousands to millions etc., of suitable compounds, for example peptides, peptoids and other oligomeric
15 compounds (cyclic or linear), and template-based smaller molecules, for example benzodiazepines, hydantoins, biaryls, carbocyclic and polycyclic compounds (eg. naphthalenes, phenothiazines, acridines, steroids etc.), carbohydrate and amino acids derivatives, dihydropyridines, benzhydryls and heterocycles (eg. triazines, indoles, thiazolidines etc.). The numbers quoted and the types of compounds listed are illustrative, but
20 not limiting. Preferred chemical libraries comprise chemical compounds of low molecular weight and potential therapeutic agents.

In a further aspect of the invention we provide a MurD enzyme modulator resulting from use of the assay method of the invention.

The activator-independent MurD enzyme may be produced using known recombinant
25 techniques for cloning and expression (cf. Sambrook, J., E. F. Fritsch, and T. Maniatis. 1989. Molecular cloning: a laboratory manual, 2nd ed. Cold Spring Harbor Laboratory, Cold Spring Harbor, N.Y.). Convenient expression systems for the MurD enzyme include T7 promoter-driven transcription of the *murD* gene in a suitable host, more conveniently *E. coli*. Examples of convenient expression vectors include those with a T7 promoter and suitable cloning sites
30 such as pET28b and pET30a (Novagen Inc. Madison WI USA). The *E. coli* host strains used for expression in such a system include those that contain the T7 RNA polymerase gene that can be induced to initiate transcription of the *murD* gene, more specifically BL21(DE3) and HMS174(DE3).

The invention will now be illustrated by reference to the following Specific Description and Figures wherein:

5 Figure 1 shows the lack of salt activation for *E. faecalis* Mur D and shows that all salts but SO_4 ions (=inhibitor) have no effect on *E. faecalis* Mur D

Figure 2 shows the salt activation of *S. aureus* MurD

10 Figure 3 shows that the *E. coli* enzyme activity is dependent on the presence of ammonium formate while the *E. faecalis* enzyme activity is independent on the presence of this salt.

Figure 4 shows the stability of *E. faecalis* Mur D in DMSO at concentrations of up to 5%.

15 Figure 5 shows a graph (Eyring Plot) with the activation enthalpies for several Mur D orthologues, including *E. faecalis*. The higher this value the more temperature dependent the reaction rate

20 Figure 6 (a) and (b) show the pH dependence of several Mur D orthologues, including *E. faecalis*

Figure 7 shows the background ATPase activity of several Mur D orthologues, including *E. faecalis*.

25 Figure 8 shows the full amino acid sequences and alignments for the Mur D enzymes from *E. coli*, *E. faecalis*, *P. aeruginosa* and *S. Aureus*.

Specific Description:

We have established that the *E. faecalis* Mur D enzyme has the following additional advantages for screening purposes.

- 5 (a) it is not affected up to a dimethylsulfoxide (DMSO) concentration of 5%. Since DMSO is commonly used in screening and IC₅₀ measurements, this contributes to the stability of the assay results.
- (b) it has a low temperature dependence. Temperature variations during assays have only a minor effect on the enzyme.
- (c) it has the broadest pH optimum (pH 7.5 – 9.5). Changes in pH due to compound addition
10 is less likely to affect activity and therefore assay results
- (d) it has a low background ATPase activity. An idle ATPase activity can create a background signal in a screen (phosphate is being produced without catalysis occurring) so a low amount of this activity is desirable to produce a sensitive assay.
- 15 We devised a screening assay for this ortholog using Malachite Green for detection of the phosphate product. The assay includes controls to check for compound interference with the signal and interference with the detection method.
 1. Compounds are preincubated with *E. faecalis* MurD for 5-30 min.
 2. The 3 substrates ATP, UDP-MurNac-L-Ala and D-Glutamate are added to initiate the
20 reaction.
 3. After 30-60 min the reaction is stopped with a Malachite Green solutionThe signal is recorded spectrophotometrically 4-10 min after the Malachite green solution had been added.
- 25 To control for compound interference with the assay signal, compounds are submitted to the same procedure but in the absence of *E. faecalis* MurD in step 1. Compounds interfering with the signal show an increased signal relative to a control where no *E. faecalis* MurD was present
- 30 To control for compound interference with the Malchite green detection method, a constant amount of phosphate (10-15 uM) is replaced for *E. faecalis* MurD in step 1. Interference is detected by a increase or decrease of the signal relative to a control where no compound was present.

100874

- 7 -

Specific assay conditions are 0.3nM *E. faecalis* MurD in 50mM Tris, 2.5mM DTT, 10mM MgCl₂, 0.01% Triton X-100, 50 uM ATP, 50uM UMA, 100uM D-Glu, pH 8.0. Typically MurD is preincubated with inhibitor for 15 min in the absence of substrates. Subsequently the reaction is initiated by adding substrates and stopped by addition of Malachite green after

5 60min. Signal is read 5 min after stopping the reaction.

100874

- 8 -

Claim:

1. The use of an activator-independent MurD enzyme in a screening assay to identify inhibitors of the enzyme, which assay comprises contacting the enzyme with a test compound
- 5 in the presence of an enzyme substrate and appropriate buffers etc. and detecting any modulation of enzyme activity by the test compound.

Figure 1

The lack of salt activation for *Enterococcus faecalis* MurD

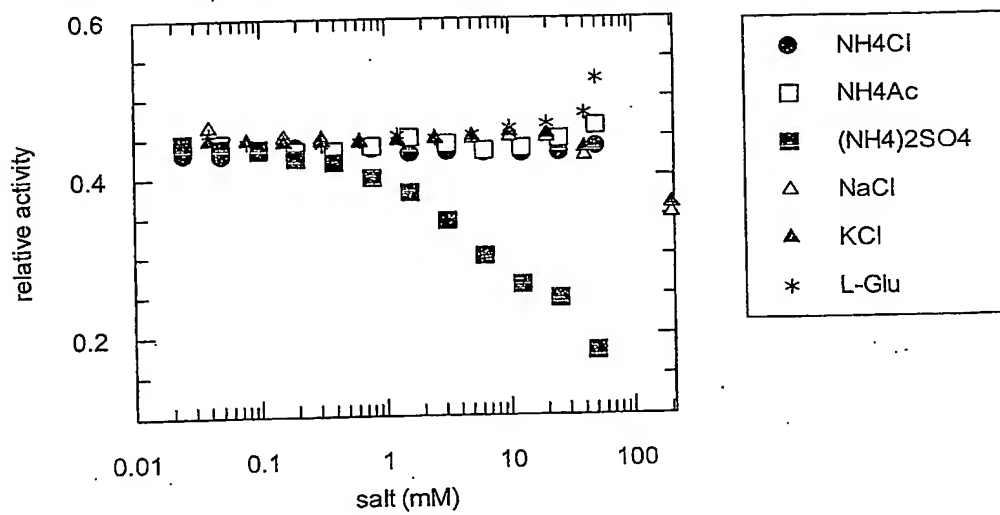


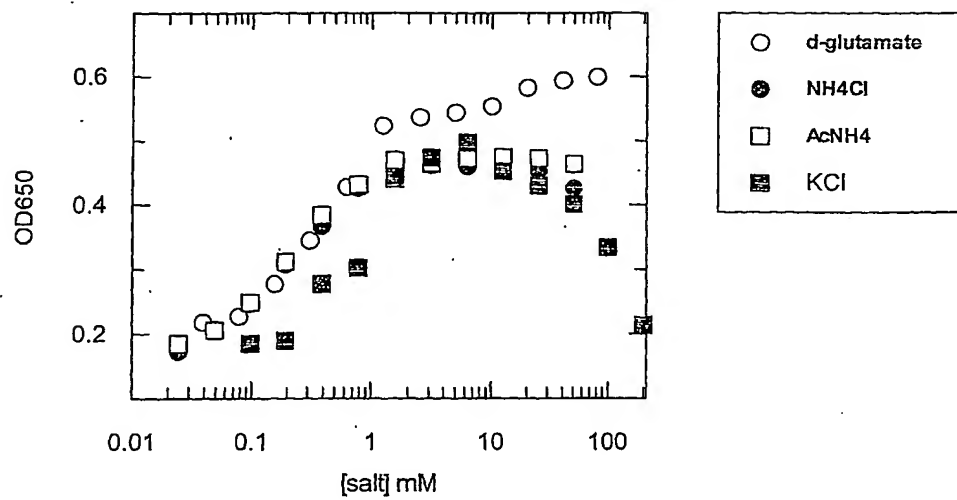
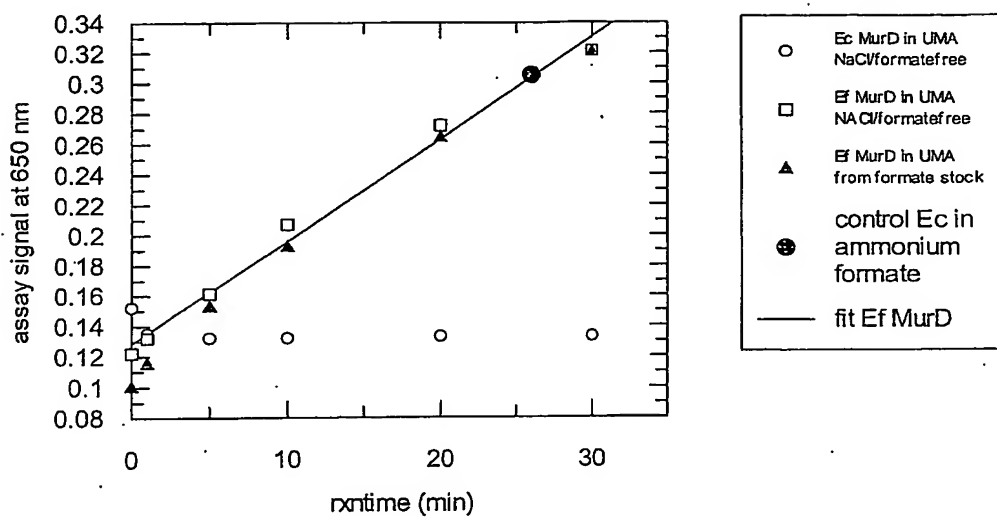
Figure 2Salt activation of *S. aureus* MurD

Figure 3



| Parameter | Value | Std. Error |
|---------------|--------|------------|
| a (intercept) | 0.1285 | 0.0054 |
| b (gradient) | 0.0067 | 0.0004 |

100874

- 4/8 -

Figure 4

DMSO dependence

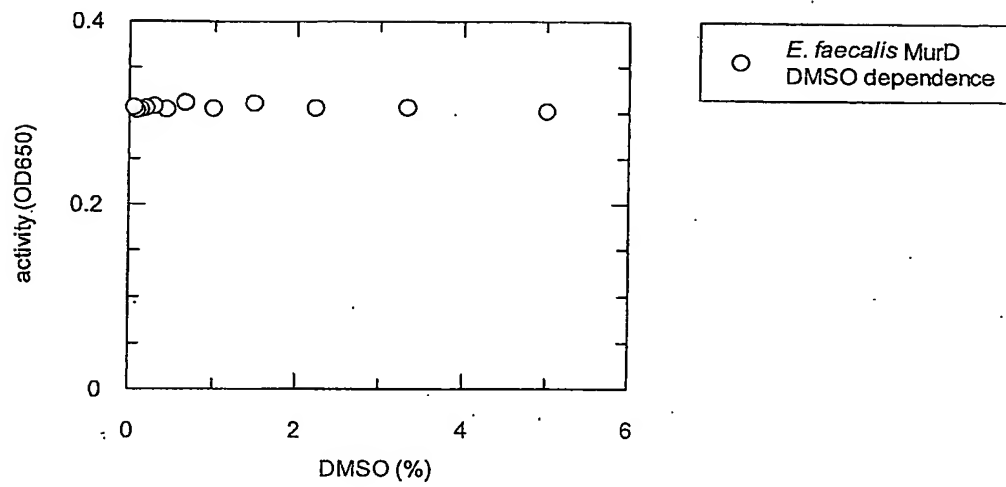


Figure 5**Temperature dependence**

This Graph (Eyring Plot) shows the activation enthalpies for all for orthologs. The higher this value the more temperature the reaction rate.

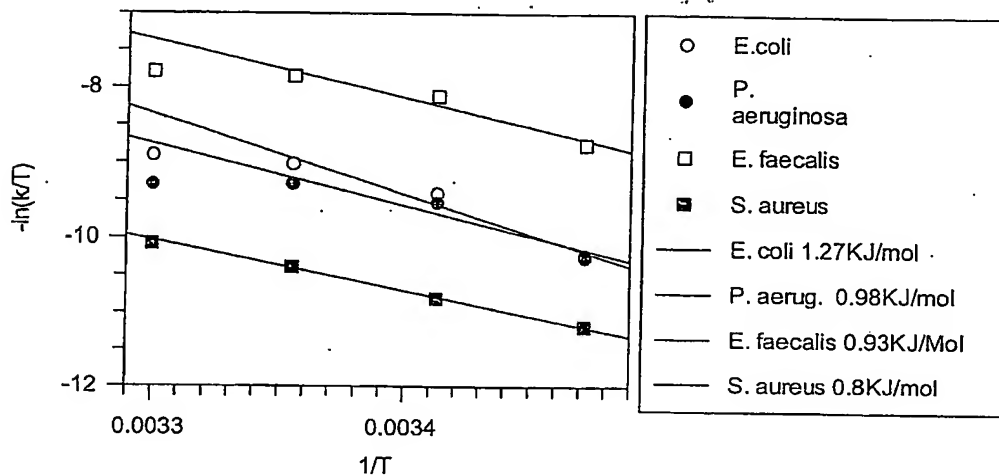
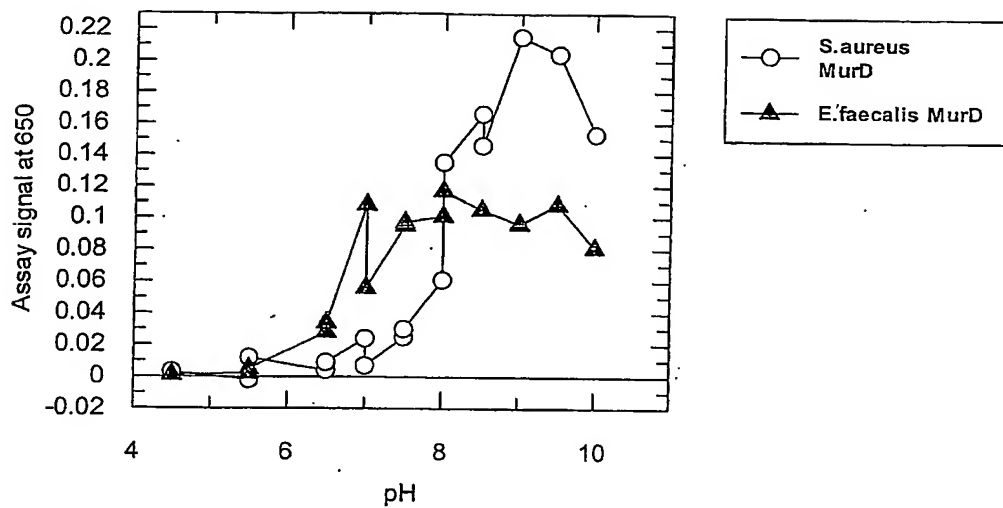


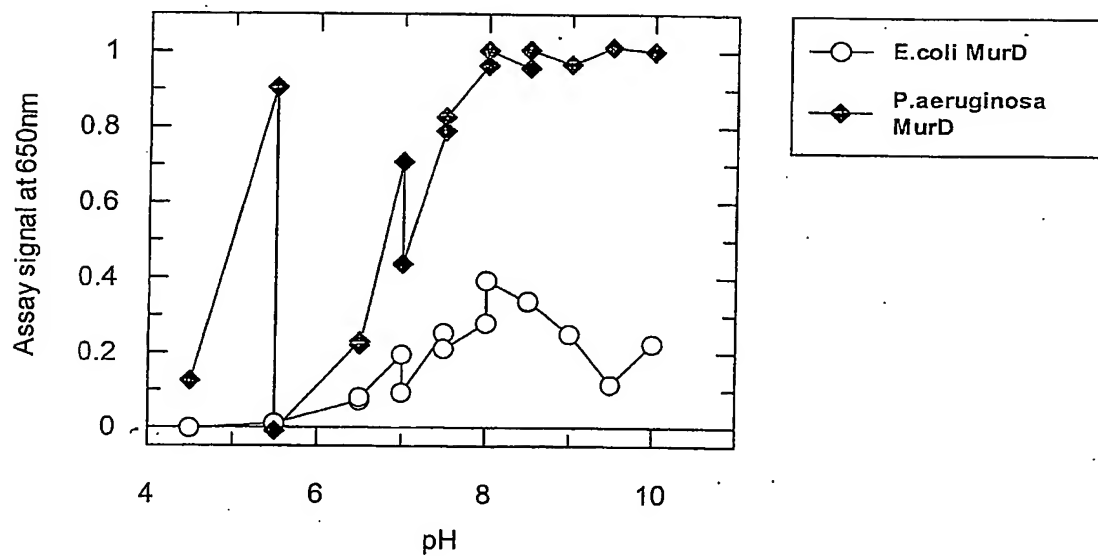
Figure 6

pH dependence

(a)



(b)

Figure 7

Background ATPase activity:

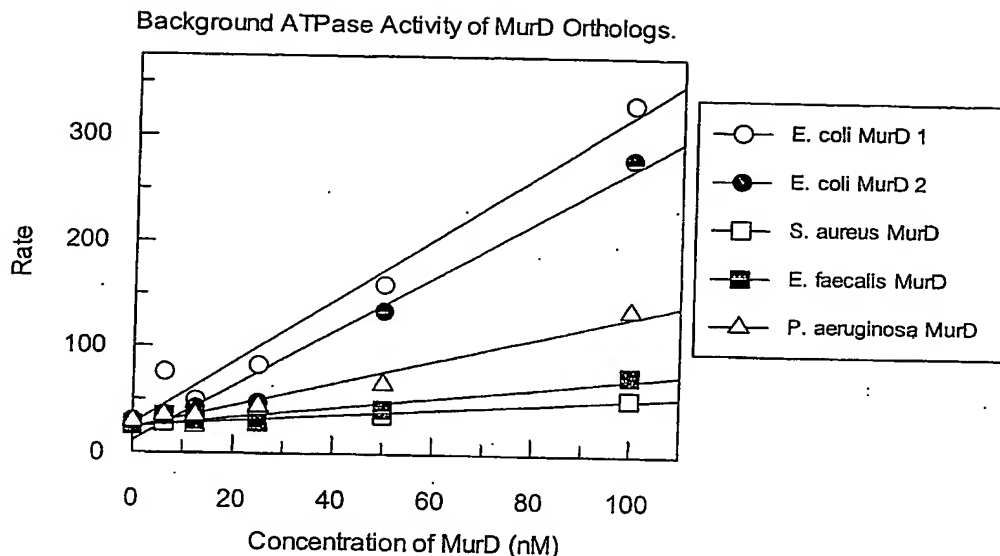


Figure 8

This alignment shows 19.2% identity and 32.6% similarity between all four MurD isozymes

```

Ecoli      ---MADYQGKNNVVIIGLGLTGLSCVDFFLARGVTPRVMDTRMTPP---GLDKLPEAVERH
Efaecalis  MKKITTYQNKVLVLGLAKSGVSAAKLLHELGAIVTVNDKQFDQNPDAQDLLTLGIRVV
Paeruginosa --MSLIASDHFRIVVGLGKSGMSLVRYLARRGLPFAVVDTRENPPELATLRAQYPQVEVR
Saureus     MLNYTGLENKNVLVGLAKSGYEAAKLLSKLGANVTVNDGKDLSQDAHAKDLESMSGISVV
           . . . . * . * . * . * . * .

Ecoli      TGSLNDEWLMAA-DLIVASPGIALAHPSLSAAADAGIEIVGDIELFCREAQAPIVAITGS
Efaecalis  TGGHPIELLDEEFELIVKNPGILYTNPLVAEALTRKIPITEVELAGQIAECPIVGITGT
Paeruginosa CGELDAEFLCSA-RELYVSPGLSLRTPALVQAAAGVVRISGDIDLFAREAKAPIVAITGS
Saureus     SGSHPLTLLDNN-PIIVKNPGIPYTVSIIIDEAVKRGKILTEVELSYLISEAPIIAVTGT
           . * . * . . . . * . . . * . . . . . . . . . . . . . .

Ecoli      NGKSTVTTLVGEMAKAAGVNV--GVGGNIGLPA--LMLLDDECELYVLELSSFQLETTSS
Efaecalis  NGKTTTTTIGLLLNADRTAGEARLAGNIGFPASTVAQEATAKDDLVMELSSFQLMGIET
Paeruginosa NAKSTVTTLVGEMAVAADKRV--AVGGNLGTPA--LDLLADDIELYVLELSSFQLETCDR
Saureus     NGKTTVTSIGDMFKKSRLTG--RLSGNIGYVASKVAQEVKPTDYLVTTELSSFQLLGIEK
           * * . * * . . . . . * * . * . * * * * * * * * * * * *

Ecoli      LQAVAATILNVTEDHMDRYPFGLQQYRAAKLRIYENAKVCVNVNADDALTMPIRG-ADERC
Efaecalis  FHPQIAVITNIFEAHLDYH-GSRKEYVAAKWAIQNMATAEDTLILNWNQVELQTLAKTTA

```

Paeruginosa LNAEVATVLNVSEDHMDRY-DGMADYHLAKHRIFRGARQVVVNRADALTRPLIA-DTVPC
 Saureus YKPHIAITNIYSAHLDYH-ENLENYQNAKKQIYKNQTEEDYLCNYHQROVIE-SEELK

* . * . * . * * *

Ecoli VSFGVN---MGDYHLNHQQGETWLRVKGEKVLNVKEMKLSGQHNYTNALALADAAGL
 Efaecalis ANVLPFSTKEAVEGAYLLDG--KLYFNEEYIMPADELGIPGSHNIENALAAICVPKLNK
 Paeruginosa WSFGLNKPDKAFGLIEEDGQKWLAFOFDKLLPVGELKIRGAHNSNALAALALGHAVGL
 Saureus AKTLYFSTQQEVDGIYIKDG--FIVYKGVRIINTEDLVLPGEHNLLENILAAVLACILAGV

* . . . * * * * *

Ecoli PRASSLKALTFTTGLPHRFEVVLEHNGVRWINDSKATNVGSTEALN--GLHVDGTLHLL
 Efaecalis SNAQIKQSLTNFSGVPHRTQFVGEVQORRFYNSKATNFLATEMALS--GFDNQKLLLLA
 Paeruginosa PFDAMLGALKAFSGLAHRQCQWVRERQGVSYDDSKATNVGAALAAIEGLGADIDGKLVLL
 Saureus PIKAIIDSLTTFSGIEHRLQYVGTNRTNKYYNSKATNTLATQFALN--SFN-QPIIWL

* * . * . * * * * * . *

Ecoli LGGDGKSAFSPRARYLNGDNVRLYCFGRDGAQLAALR--PEVAEQT--ETMEQAMRLLA
 Efaecalis GGLD-RGNSFDELVPALLG-LKAIIVLFGETKKKLAEAAKPNITILFAENVQTAVTIAF
 Paeruginosa AGGDGKGADFHDLEPVARFCRAVVLLGRDAGLIAQAL--GNAVPLVRVATLDEAVRQAA
 Saureus GGLD-RGNEFDELIPYMEN-VRAMVVFGQTKAKFAKLNS-QGKSVIEANNVEDAVDKVQ

* * . * * . * * * * . *

Ecoli PRVQPGDMVLLSPACASLDQFKNFEQRGNEFARLAKELG----
 Efaecalis DYSEKDDTILLSPACASWDQYPNFEVRGEAFMOAVQQLKESEM
 Paeruginosa ELAREGDAVLLSPACASLDMFKNFEERGRLFAKAVEELA----
 Saureus DIIEPNDVVLLSPACASWDQYSTFEERGEKFTIERFRAHLPSY-

* * * * * * * * * *

PCT Application
GB0304592



**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ BLACK BORDERS
- ☐ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES
- ☐ FADED TEXT OR DRAWING
- ☐ BLURRED OR ILLEGIBLE TEXT OR DRAWING
- ☐ SKEWED/SLANTED IMAGES
- ☐ COLOR OR BLACK AND WHITE PHOTOGRAPHS
- ☐ GRAY SCALE DOCUMENTS
- ☒ LINES OR MARKS ON ORIGINAL DOCUMENT
- ☐ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY
- ☐ OTHER: _____

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.